BLACK

7.16 mm.

SIZE : 148x105 mm.

Panasonic

Operating Instructions
(Household use) Rechargeable Oral Irrigator

Model No. EW-DJ67
EW-DJ66

Series 600

Panasonic Corporation
<https://www.panasonic.com>
© Panasonic Corporation 2024

© EN, DE, FR, IT, SP, NL, PT, SV, DA, NB, FI, PL, CS, SK, HU, RO, TR

Printed in Thailand
SEWDJ66303U03 A0424-0

Please check here for Instructions
and Troubleshooting.

PANASONIC : SEWDJ66303U03_COVER_OL_PROOF2_NUK_28/2/2024
PANASONIC : SEWDJ66303U03_COVER_ORIGINAL_PROOF1_NUK_28/2/2024

ILLUSTRATOR OK

1
 2
 3
 4
 5
 6

A blank coordinate plane with a horizontal x-axis and a vertical y-axis intersecting at the origin. The axes are represented by thin black lines.


(i)

	C	(iv)	D	E
EW-DJ67				
EW-DJ66				

A diagram of a corded electric kettle. Label (ii) points to the power cord, and label (iii) points to the three-prong electrical plug.

1
 2
 3
 4
 5
 6

ILLUSTRATOR OK

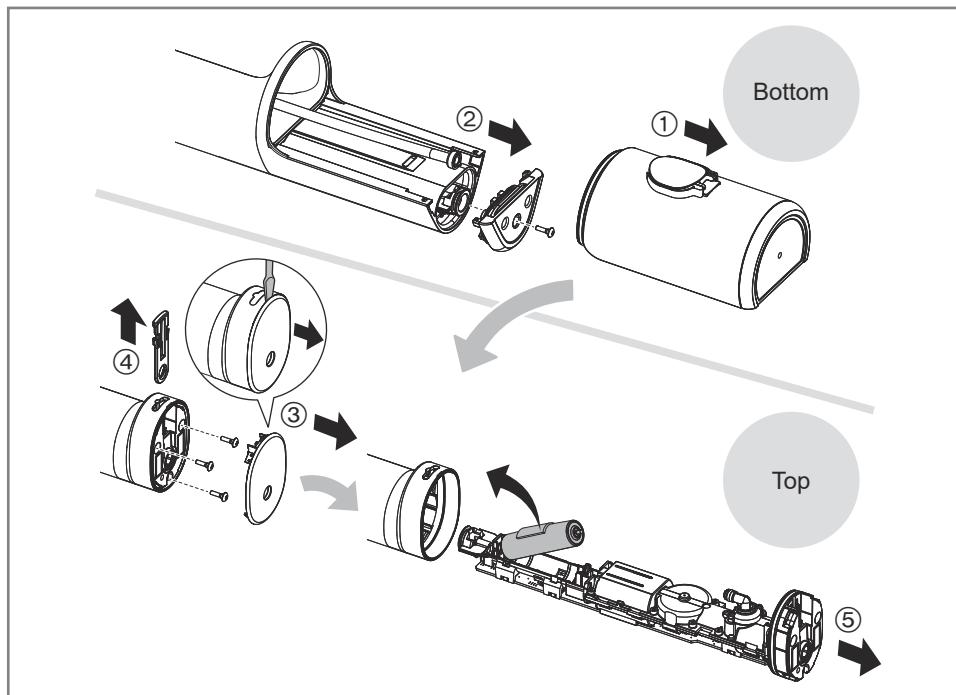English	5	Português	59	Česky	109
Deutsch	13	Svenska	68	Slovensky	117
Français	23	Dansk	76	Magyar	125
Italiano	32	Norsk	84	Română	134
Español	41	Suomi	92	Türkçe	143
Nederlands	50	Polski	100		

Panasonic

English

Operating Instructions (Household use) Rechargeable Oral Irrigator

Model No. EW-DJ67
EW-DJ66

Contents

Safety precautions.....	6	Battery life	10
Important information	9	Disposal of the appliance	10
Parts identification	10	Specifications	11
How to use	10		

Thank you for purchasing this Panasonic product.

Before operating this unit, please read these instructions completely and save them for future use.

Please read the code on the cover for the detailed operating instructions.

WARNING

- Do not immerse the charger in water or wash it with water.
Doing so may cause electric shock or fire due to a short circuit.
- Always unplug the power plug from a household outlet when cleaning the charger.
Failure to do so may cause electric shock or injury.
- The supply cord cannot be replaced. If the cord is damaged, the charger should be scrapped.
Failure to do so may cause an accident or injury.
- Oral irrigators can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge under supervision.
- Oral irrigators shall be used only with cold or warm water, or specific solutions as defined in the manufacturer's instructions.
- Children shall not play with the appliance.
- Do not use anything other than the supplied charger for any purpose. Also, do not use any other product with the supplied charger.
Doing so may cause burns or fire due to a short circuit.
- The following symbol indicates that a specific detachable power supply unit is required for connecting the electrical appliance to the supply mains. The type reference of power supply unit is marked near the symbol.

Safety precautions

Make sure to follow these instructions.

To prevent accidents, injuries or damage to property, please follow the instructions below.

- The following chart indicates the degree of damage caused by improper operation.

DANGER

Indicates a potential hazard that will result in serious injury or death.

WARNING

Indicates a potential hazard that could result in serious injury or death.

CAUTION

Indicates a potential hazard that could result in minor injury or property damage.

WARNING

► Preventing accidents

- Do not store within the reach of children or infants. Do not let them use it.
- Doing so may cause an accident or injury due to accidental ingestion of removable parts.

 WARNING**► Power supply**

- Do not connect or disconnect the power plug to a household outlet with a wet hand.
 - Doing so may cause electric shock or injury.
- Do not place the charger over or near water filled sink or bathtub.
- Do not use when the power cord or the power plug is damaged or when the fitting into the household outlet is loose.
- Do not damage, modify, forcefully bend, pull, twist, or bundle the power cord or the power plug. Also, do not place anything heavy on or pinch the power cord.
 - Doing so may cause electric shock or fire due to a short circuit.
- Do not use in a way exceeding the rating of the household outlet or the wiring.
 - Exceeding the rating by connecting too many power plugs to one household outlet may cause fire due to overheating.
- Always ensure the appliance is operated on an electric power source matched to the rated voltage indicated on the charger.
- Fully insert the power plug.
 - Failure to do so may cause fire or electric shock.

► In case of an abnormality or malfunction

- Immediately stop using and unplug the power plug when using the charger.
 - Failure to do so may cause fire, electric shock or injury.
 - <Abnormality or malfunction cases>
 - The main unit or charger is deformed or abnormally hot.
 - The main unit or charger smells of burning.
 - There is abnormal sound during use or charging of the main unit or charger.
- Immediately request inspection or repair to an authorized service centre.

► This product

- This product has a built-in rechargeable battery. Do not throw into fire, apply heat, or charge, use, or leave in a high temperature environment.
 - Doing so may cause overheating, ignition, or explosion of the built-in rechargeable battery.
- Do not modify or repair.
 - Doing so may cause fire, electric shock, or injury.
- Contact an authorized service centre for repair (battery change, etc.).
- Never disassemble except when disposing of the product.
 - Doing so may cause fire, electric shock, or injury.

► Cleaning

- Regularly clean the power plug to prevent dust from accumulating (once every half a year).
 - Failure to do so may cause fire due to insulation failure caused by humidity.
- Disconnect the power plug and wipe with a dry cloth.

**CAUTION****► Preventing damage to teeth and gums**

- People who are unable to operate should not use this appliance.
- People with weak oral senses should not use this appliance.
 - Doing so may cause injury or damage to teeth, gums and tongue.
- People who may have periodontal disease, have teeth being treated, or concerned about symptoms within their mouth should consult a dentist before use.
- Observe the following precautions when using the tongue care nozzle.
- Do not allow children to use it.
- Do not brush forcefully.
- Do not insert too deep (approx. 3/4 of the visible range).
- If tongue inflammation occurs, stop using the device and consult a doctor if the inflammation persists after 3 days.
 - Failure to do so may cause injury or damage to teeth, gums and tongue.

► This product

- Never use for a purpose other than oral cleaning.
 - Doing so may cause an accident or injury when the stream of water is pointed toward the eyes, nose, ears, tongue, or throat.
- Do not fill the water tank with water warmer than 40 °C.
 - Doing so may cause burns.

- Do not share nozzle with your family or other people.
 - Doing so may cause infection or inflammation.
- Do not place a metallic object such as a coin or a clip on top of the charger.
 - Doing so may cause burns due to excessive heat.

► Note the following precautions

- Do not allow metal objects or trash to adhere to the power plug.
 - Doing so may cause electric shock or fire due to a short circuit.
- Do not drop or subject to shock.
 - Doing so may cause injury.
- Do not wrap the power cord around the charger when storing.
 - Doing so may cause the wire in the power cord to break with the load, and may cause fire due to a short circuit.
- Disconnect the power plug from the household outlet when not charging.
 - Failure to do so may cause electric shock or fire due to electrical leakage resulting from insulation deterioration.
- Disconnect the power plug by holding onto it instead of the power cord.
 - Failure to do so may cause electric shock or fire due to a short circuit.
- After use, clean the main unit and the nozzle and dry them.
 - Failure to do so may cause damage to your health due to propagation of mold and bacteria.

Handling of the removed battery when disposing

DANGER

- The rechargeable battery is exclusively for use with this appliance. Do not use the battery with other products.
- Do not charge the battery after it has been removed from the product.
- Do not throw into fire or apply heat.
- Do not hit, disassemble, modify, or puncture with a nail.
- Do not let the positive and negative terminals of the battery get in contact with each other through metallic objects.
- Do not carry or store the battery together with metallic jewellery such as necklaces and hairpins.
- Do not use or leave the battery where it will be exposed to high temperatures, such as under direct sunlight or near other sources of heat.
- Never peel off the tube.
 - Doing so may cause overheating, ignition, or explosion.

WARNING

- After removing the rechargeable battery, do not keep it within reach of children and infants.
 - The battery will harm the body if accidentally swallowed.
- If this happens, consult a physician immediately.

- If the battery fluid leaks out and comes in contact with the eyes, do not rub the eyes, and thoroughly rinse with clean water, such as tap water.
 - Failure to do so may cause eye injury.
- Consult a physician immediately.

CAUTION

- If the battery fluid leaks out and comes in contact with the skin or clothes, rinse off with clean water, such as tap water.
 - Failure to do so may cause inflammation.

Important information

- When you first begin using the oral irrigator, you may experience slight bleeding of the gums even though your gums may be healthy. This is because your gums are receiving stimulation by the oral irrigator for the first time. The bleeding should stop after one or two weeks of use. Bleeding that persists for more than two weeks may be indication of an abnormal condition of your gums, such as periodontal disease. In this case, stop using the oral irrigator and consult your dentist.
- Do not put salt, mouthwash or other chemicals into the water tank. Doing so may cause malfunction.
- Do not dry the appliance with a dryer or a fan heater. Doing so may cause malfunction or deformation of the parts.

- Store the appliance in a place with low humidity after use. Leaving the appliance in a bathroom may cause malfunction.
- Charge the appliance correctly according to these Operating Instructions.
- This appliance contains a battery that is only replaceable by skilled persons. Please contact an authorized service centre for the details of repair.

Parts identification (page 2)

- A Main unit**
 - (i) Water tank
- B Charger (RE8-59)**
 - (ii) Power cord
 - (iii) Power plug
- C Ultrasonic nozzle**
 - (iv) Identification ring
- D Orthodontic nozzle**
- E Tongue care nozzle**

- The included nozzles vary by model. (See table.)

How to use

Please refer to the **Quick Start Guide** for how to use the appliance.

Battery life

Approximately 3 years. The battery may have reached its operating life if the usage frequency decreases drastically, even after a full charge.
(The operating life of the battery may differ depending on usage or storage conditions.)

Disposal of the appliance

Remove the built-in rechargeable battery before disposing of the appliance.

- The battery is to be disposed of safely.
- Please make sure that the battery is disposed of at an officially designated location if there is one.

How to remove the battery (page 3)

This figure must only be used when disposing of the appliance, and must not be used to repair it. If you dismantle the appliance yourself, it will no longer be waterproof, which may cause it to malfunction.

- Remove the main unit from the charger, press the power switch to turn on the power and then keep the power on until the battery is completely discharged.
- Perform steps ① to ⑤ to disassemble the main unit and remove the battery with a screwdriver.
- Insulate the + and - terminals by applying tape to them. (To prevent short circuit of terminals).

For environmental protection and recycling of materials

This appliance contains a Lithium-ion battery.
Please make sure that the battery is disposed of at an officially assigned location, if there is one in your country.

Specifications

Power source	See the name plate on the product (Automatic voltage conversion)
Charging time	Approx. 17 hours
Operating time	LEVEL 5: Approx. 10 minutes (at 20 °C when fully charged)
Tank volume	Approx. 200 mL (For approx. 60 seconds of operation when the water pressure level is set to 5 when using the ultrasonic nozzle)
Dimensions	Main unit: 219 (H) x 63 (W) x 80 (D) mm (Not including nozzle) Charger: 37 (H) x 86 (W) x 101 (D) mm
Mass	Main unit: Approx. 275 g (Including the ultrasonic nozzle) Charger: Approx. 200 g
Airborne Acoustical Noise	60 (dB (A) re 1 pW)

This product is intended for household use only.

Authorized Representative in EU:

Panasonic Marketing Europe GmbH
Panasonic Testing Centre
Winsbergring 15, 22525 Hamburg, Germany

Importer:

Panasonic Marketing Europe GmbH
Hagenauer Strasse 43, 65203 Wiesbaden, Germany

Manufactured by:

Panasonic Corporation
1006, Oaza Kadoma, Kadoma City, Osaka, Japan

Disposal of Old Equipment and Batteries
Only for European Union and countries with recycling systems

These symbols on the products, packaging, and/or accompanying documents mean that used electrical and electronic products and batteries must not be mixed with general household waste.

For proper treatment, recovery and recycling of old products and batteries, please take them to applicable collection points in accordance with your national legislation.

By disposing of them correctly, you will help to save valuable resources and prevent any potential negative effects on human health and the environment. For more information about collection and recycling, please contact your local authority.

Penalties may be applicable for incorrect disposal of this waste, in accordance with national legislation.

Note for the battery symbol (bottom symbol)

This symbol might be used in combination with a chemical symbol. In this case it complies with the requirement set by the Directive for the chemical involved.

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Ltd.**
AB Yetkili Temsilcisi: