



SKG

Neck Massager
G7Pro-Fold

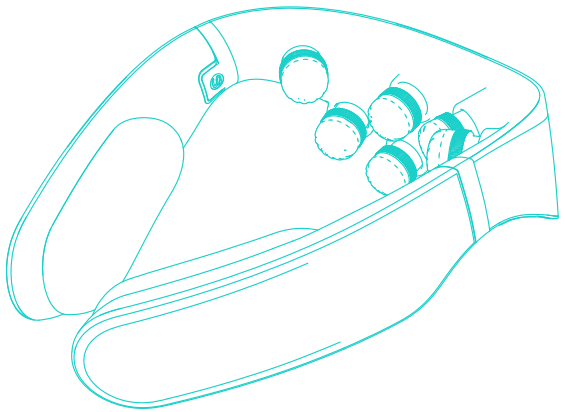

USER MANUAL


Enjoy the power-on experience

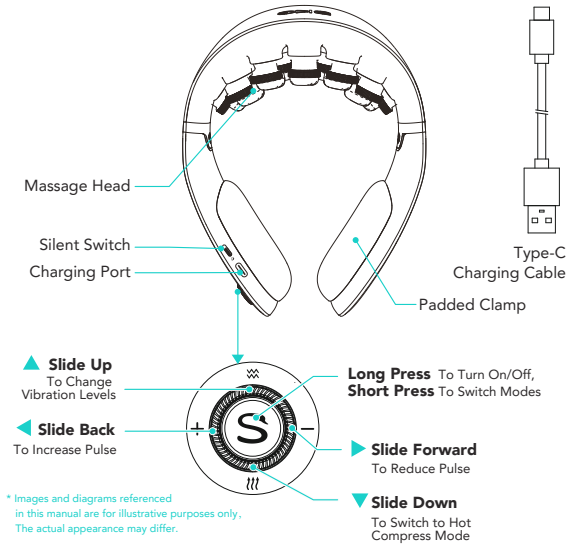

Charging/Power Instructions		
Light Color	Light Status	Description
Blue Light	Always on	Battery Level Normal
Orange Light	Always on	Low Battery Level
Orange Light	Breathing	Charging
Blue Light	Breathing	Battery Fully Charged

Basic Product Parameters			
Model	G7Pro-Fold	Rated Input Voltage	5V ^{***}
Battery Capacity	2400mAh/8.88Wh	Rated Input Power	10W


Ciesz się funkcjonalnością naszych produktów

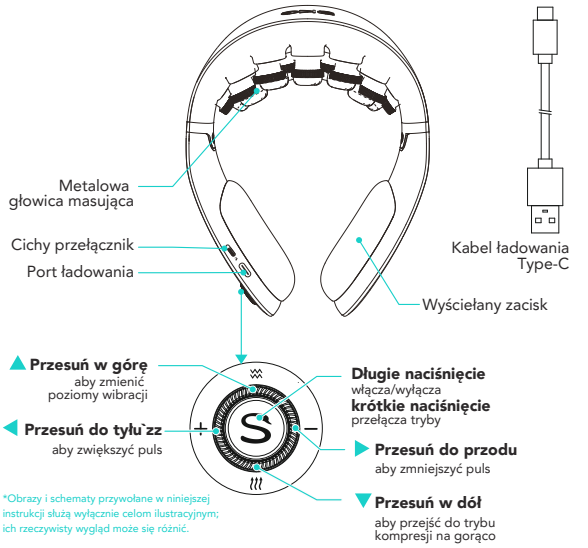

Instrukcja ładowania/zasilania		
Kolor światła	Status światła	Opis
Niebieskie światło	Zawsze włączone	Normalny poziom naładowania baterii
Pomarańczowe światło	Zawsze włączone	Niski poziom naładowania baterii
Pomarańczowe światło	Powolne miganie	Ładowanie
Niebieskie światło	Powolne miganie	Bateria w pełni naładowana

Podstawowe parametry produktu			
Model	G7Pro-Fold	Znamiennowe napięcie ładowania	5V ^{***}
Pojemność baterii	2400mAh/8.88Wh	Znamiennowa moc wyjściowa	10W


Genießen Sie das Erlebnis beim Einschalten

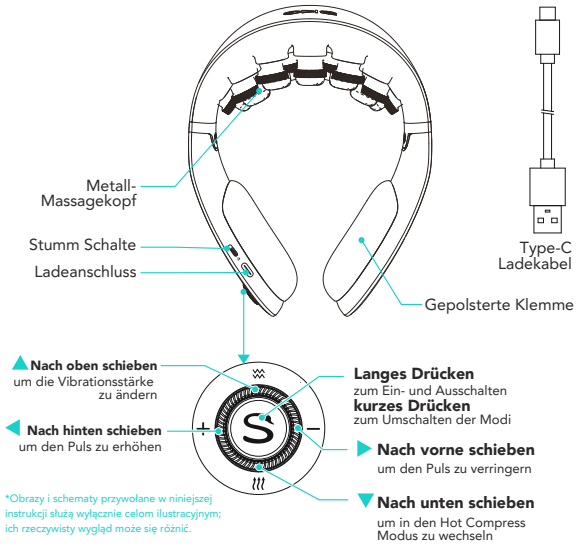

Anweisungen zum Aufladen/Stromversorgung		
Lichtfarbe	Licht-Status	Beschreibung
Voyant bleu	Immer an	Normaler Batteriestand
Oranges Licht	Immer an	Niedriger Batteriestand
Oranges Licht	Blinken	Ladet auf
Voyant bleu	Blinken	Voll aufgeladen

Grundlegende Produktparameter			
Modell	G7Pro-Fold	Nenn-Ladespannung	5V ^{***}
Batteriekapazität	2400mAh/8.88Wh	Nenn-Eingangsleistung	10W

Profitez de l'expérience de mise en marche

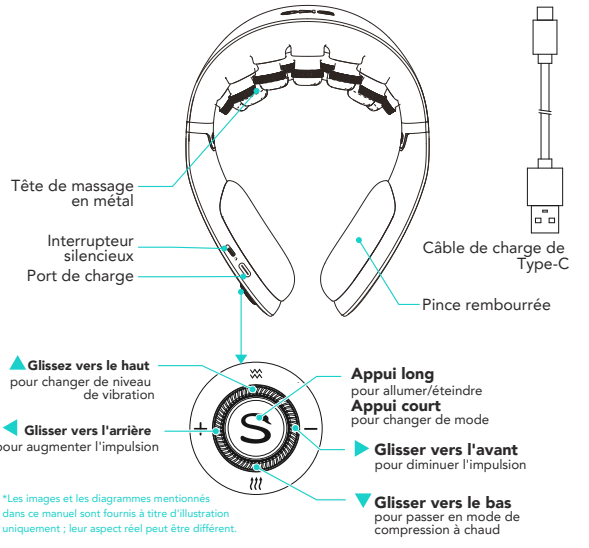

Instructions de charge et d'alimentation		
Couleur de la lumière	Statut de lumière	Description
Voyant bleu	Toujours activé	Niveau de batterie normal
Lumi ère orange	Toujours activé	Niveau de batterie faible
Lumi ère orange	Respiration	Chargement
Voyant bleu	Respiration	La batterie est complètement chargée

Grundlegende Produktparameter			
Modèle	G7Pro-Fold	Tension de charge nominale	5V ^{***}
Capacité de la batterie	2400mAh/8.88Wh	Puissance d'entrée nominale	10W


Disfruta de la mejor experiencia al encenderlo



Instrucciones de carga/encendido		
Color de luz	Estado de la luz	Descripción
Luz azul	Siempre encendido	Nivel de batería normal
Luz naranja	Siempre encendido	Nivel bajo de batería
Luz naranja	Respiración	Cargando
Luz azul	Respiración	Batería completamente cargada

Características básicas del producto			
Modelo	G7Pro-Fold	Voltaje de carga nominal	5V ^{***}
Capacidad de la batería	2400mAh/8.88Wh	Potencia de entrada nominal	10W

Godetevi un'esperienza che vi dà la carica

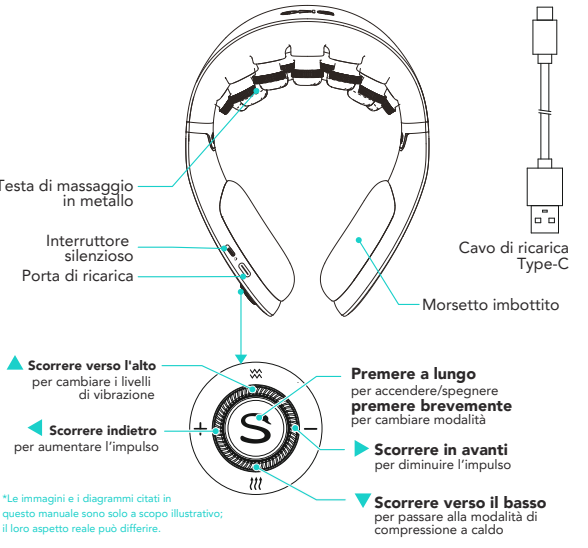

Istruzioni per la ricarica e l'alimentazione		
Colore spia	Stato spia	Descrizione
Luce blu	Fissa	Livello batteria normale
Luce arancione	Fissa	Livello basso della batteria
Luce arancione	Lampeggiante	Ricarica in corso
Luce blu	Lampeggiante	Batteria completamente carica

Parametri di base del prodotto			
Modello	G7Pro-Fold	Tensione nominale di rete	5V ^{***}
Capacità della batteria	2400mAh/8.88Wh	Potenza nominale in ingresso	10W


SKG Health Technologies Co., Ltd.
Address: 23A Floor, Building 3, Zhongke R&D Park, No. 009,
Gaixin South 1st Road, High-tech Zone Community, Yuehai Street,
Nanshan District, Shenzhen, Guangdong, China.
Website : <https://global.skg.com/>



Made in China

Wearing Method

- Please wear close to the skin
- Moisten skin with wet wipe before use

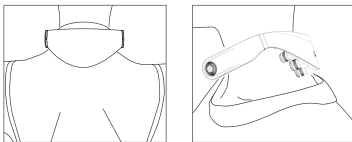

Daily maintenance of the product

- Avoid direct sunlight
- Avoid high temperatures and humidity
- Avoid being near open flames
- Keep away from hard objects, avoid knocking, and do not scratch.
- The device's exterior may be cleaned using a moist or soft cloth.
- It is recommended to charge once a month if not in use for an extended period of time and keep the battery fully charged when using

Troubleshooting Instructions

If the product is not working properly, please check for the following conditions. If our suggestions do not solve the problem, please return the product to the dealer for repair.

Condition	Possible Cause	Solution
Unable to start	Low battery power	Charge the product
	The product is not turned on	Long-press the power button to turn on
Unable to charge or no charging light displayed	The product is still connected to the household power socket	Remove the charging cable from the device
	Charging cable is not properly inserted into the power socket	Replug the charging cable
	Wrong charging cable used	Use a Type-C connector cable

Product Disposal Precautions

Please support the recycling of rechargeable batteries (lithium-ion batteries)	
Danger	This product's battery is a dedicated rechargeable battery that should not be used in other products. Do not charge the battery once it has been removed. - Do not throw into the fire or expose to excessive heat. - Do not puncture with nails, subject to heavy impact, or disassemble - Do not directly connect the battery's positive and negative terminals with metal - Do not store it with necklaces, hairpins, etc. - Do not charge, use, or store in high-temperature locations such as near a stove or under direct sunlight - All of the abovementioned scenarios can result in overheating, fire, or bursting.
Warning	Place the removed rechargeable battery out of the reach of infants to prevent accidental ingestion - If swallowed, seek medical attention immediately. - If battery leakage occurs, do not touch the liquid with bare hands. - Follow the safety instructions below: - There is the risk of blindness if battery fluid contacts the eyes. Do not rub the eye with your hands. Rinse with clean water immediately and seek medical attention. If battery liquid gets on the body or clothing, it may cause skin inflammation or injury. Rinse with clean water, and if discomfort persists, seek medical attention immediately. - Disassembly is only permitted when the product is being discarded. Failure to comply may result in fire, electric shock, injury, or other accidents. Please contact the retailer where you purchased the product for repairs (except battery replacement).

The lithium-ion battery used in this product is a valuable recyclable resource. Please remove the lithium-ion battery prior to disposal and hand over it to a professional organization capable of recycling and processing rechargeable batteries.

- The battery must be removed from the device before discarding the appliance.
- The device must be turned off before removing the battery.
- The battery should be disposed of safely.

Removal Steps:

- Open the shell using a tool.
- Unscrew the screws at the circuit board.
- Disconnect the battery terminal plug.
- Remove the battery.

Precautions

- Cautious**
- To mitigate its impacts and possible harm, do not use the product with the following medical electronic devices:
- Implanted medical electronic devices such as pacemakers.
 - Life-sustaining medical electronic devices such as artificial hearts and lungs.
 - Medical electronic devices such as electrocardiograms.
 - Use is prohibited in flammable and explosive areas such as gas stations, high-dust workshops, and gas stations.
 - This product is not equipment and does not diagnose or treat neck disorders. Patients with neck ailments should consult healthcare professionals to avoid delaying treatment.
 - Do not use this product while driving, climbing, or engaging in other dangerous activities, as it may lead to accidents.
 - Do not use near the heart, above the neck, head, mouth, genital areas, or on skin disorders, and avoid massaging both soles of the feet simultaneously.
 - The hot compress period for this product must not exceed 15 minutes (please confirm the time), as there is a risk of low-temperature burns. Do not use while asleep or in an unconscious state, and users who are sensitive to heat should use a lower heat setting and strictly control the use time.
 - Individuals with circulatory disorders such as cerebral vascular diseases, inflamed or swollen neck skin, suppurative, or neck-related ailments should not use this product, as it may worsen symptoms.
 - Keep this product out of reach of infants, children, or individuals who are entirely or partially unconscious. Do not let them use it, as it may lead to accidents.
 - Individuals over 70 years old are prohibited from using this product. Please select a lower setting and strictly control the usage duration since those over 65 have more delicate skin and less sensitive nerve endings.
 - Gradually increase the massage setting according to personal comfort; otherwise, the discomfort from a higher setting might cause muscle strain.

- Warning**
- This is a personal care massage product and cannot be used to substitute medicines or medical devices.
 - Children, pregnant women, and women on their menstrual period are prohibited from using this product.
 - Contains a lithium battery. Do not place it in a high-temperature environment.
 - Use this product under the supervision of a doctor if you have a sensitive skin, facial allergic dermatitis, or facial areas containing cosmetic materials.
 - Users allergic to metal or heat must take precautions.
 - This product has a heating function; please be cautious if you are insensitive to heat.

- If discomfort occurs during use, please stop using immediately.
- Do not disassemble this product yourself.
- Do not attempt to replace the battery yourself.
- If you have a neck injury, consult healthcare professionals before using this product.
- The device has a heating surface, so those who are insensitive to heat must take caution while using it.
- Use within the normal operating environment: Temperature: 0°C to 40°C; Relative Humidity: <80%RH, no condensation; Atmospheric Pressure: 50Kpa to 106Kpa.
- Please choose normal transportation and storage conditions: Temperature Range: -20°C to +55°C; Relative Humidity Range: ≤95%RH, no condensation; Atmospheric Pressure: 50Kpa to 106Kpa.
- This product contains a lithium battery. Do not disassemble, modify, immerse in water, or expose to open flame, as it may result in damage or explosion.
- Do not press on the heart when using, as this may cause an accident or physical discomfort.
- Do not use this product with medical devices or with more than two of this product simultaneously, as it may cause nausea or physical discomfort.
- Do not use for purposes other than those specified in the instruction manual, as it may cause accidents.
- Do not disassemble, repair, or modify this product yourself, as it may lead to a malfunction or an electric shock accident.
- Do not use to massage parts of the body other than the neck.
- If a patient simultaneously uses high-frequency surgical equipment and this massager, it may cause burns at the massager's electrode location and may also damage this massager.
- If this product is used near shortwave or microwave therapy equipment (within 1 meter), the output of this massager may be unstable.
- Electrode pads should avoid contact with wounds and scars during massage.
- While massaging, the electrode pads must not be placed near the heart area of the chest, as this could increase the risk of cardiac fibrillation.

- In the following situations, please consult with a physician before use, otherwise accidents or discomfort may occur:**
- Patients currently under medical treatment and experiencing abnormal physical symptoms.
 - Patients with infectious diseases.
 - Patients allergic to stainless steel electrode pads.
 - Patients with heart or brain nerve abnormalities.
 - Patients with sensory impairments or skin abnormalities.
 - Patients with a body temperature above 38 degrees Celsius (during a fever).
 - Patients taking medications.
 - Patients with abnormal blood pressure.
 - Patients who are unable to express their own will freely.
 - Patients with perception disorders due to severe peripheral circulation problems caused by conditions such as diabetes.

- The following situations prohibit use:**
- Patients with malignant tumors.
 - Pregnant women and women in childbirth.
 - Patients with acute suppurative inflammation.
 - Individuals with severe cognitive impairments.
 - Persons with spastic paralysis.
 - Individuals with abnormal nervous stress.
 - Individuals in a spasm phase following a cerebrovascular accident.
 - Individuals with unclear consciousness, and limb and joint contracture deformities.
 - Patients with a tendency to bleed.
 - Users of cardiac pacemakers.
 - Eye care
 - Patients with metal implants.
 - Individuals with clotting disorders.
 - Near the human heart.

- Note**
- Please read the instruction manual thoroughly before using this product.
 - If you feel physical discomfort due to abnormalities in this product, please stop using it immediately. Consult a doctor if you are experiencing unusual sensations in your body or skin, and follow medical advice.
 - Do not massage with electronic devices, as it may cause malfunctions (such as a disruption in the electronic watch's timing).
 - If a fault occurs and the device cannot operate, immediately disconnect the power and remove it; otherwise, it may cause overheating or short circuits, resulting in an accident. Contact customer service for after-sales consultation.
 - Do not plug or unplug the charging cable with wet hands, as it may cause electric shock or accidents.
 - Do not use while sleeping, as it may cause the device to malfunction and the electrode pads to move to other parts of the body, causing discomfort.
 - Please wear this product correctly for the cervical spine area, placing it symmetrically on both sides to avoid accidents or discomfort.
 - When disconnecting the charging cable from the device, hold the plug rather than the cable, since pulling may cause an electric shock or damage the product.
 - Charge the product according to the instruction manual the first time you use it or if the battery has been drained from long-term storage to ensure it works properly.
 - It is normal for the product's surface to get heated after extended use.
 - When the product is working, it will moderately knead the neck area; do not massage for an extended period of time. If discomfort is felt, take a break before using again.
 - Do not forcefully fold the product when using or storing it to prevent damaging its components.
 - During the folding process of the product, avoid hair or skin contact to prevent pinching.

Please do not use or store this product in the following locations:

- In places exposed to direct sunlight
- In strong electromagnetic fields
- In hot and humid locations
- In locations susceptible to vibration and shock
- In places with water
- Near fireworks
- In places with excessive dust

- Recommendation**
- Prior to use, please ensure the main host buttons are working properly.
 - If using the product for the first time or after the battery has been depleted from long-term storage, charge the product according to the instruction manual to ensure it works properly.
 - Please dispose of this product in line with applicable city environmental protection standards to prevent environmental pollution.
 - Recommended for use by individuals aged 16 and above.

This device complies with part 15 of the FCC Rules. Operation is subject to the condition that this device does not cause harmful interference (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.

RSS-Gen Issue 4 December 2014 * & "CNR-Gen 4e Décembre 2014:

- English: This device complies with Industry Canada licence-exempt RSS standard(s). Operation is subject to the following two conditions: (1) This device may not cause interference, and (2) This device must accept any interference, including interference that may cause undesired operation of the device.
- French: Le présent appareil est conforme aux CNR d'Industrie Canada applicables aux appareils radio exempts de licence. L'exploitation est autorisée aux deux conditions suivantes: (1) l'appareil ne doit pas produire de brouillage, et (2) l'utilisateur de l'appareil doit accepter tout brouillage radioélectrique subi, même si le brouillage est susceptible d'empêcher le fonctionnement.

CE This device was tested for typical body-worn operations. To comply with RF exposure requirements, a minimum separation distance of 0cm must be maintained between the user's body and the handset, including the antenna. Third-party belt-clips, holsters, and similar accessories used by this device should not contain any metallic components. Body-worn accessories that do not meet these requirements may not comply with RF exposure requirements and should be avoided. Use only the supplied or an approved antenna. This device is in compliance with the essential requirements and other relevant provisions of Directive 2014/53/EU. All essential radio test suites have been carried out.

1. CAUTION - RISK OF EXPLOSION IF BATTERY IS REPLACED BY AN INCORRECT TYPE. DISPOSE OF USED BATTERIES ACCORDING TO THE INSTRUCTIONS

2. The device complies with RF specifications when the device used at 0mm from your body EU Compliance Statement:

SKG Health technology Co., Ltd. Here by declares that this device is in compliance with essential requirements and other relevant provisions of the Directive 2014/53/EU. A copy of the EU Declaration of Conformity is available online at <https://www.skg.com/pages/skg-declaration-of-conformity>.

Old electrical appliances must not be disposed of together with the residual waste, but have to be disposed of separately. The disposal at the commune collecting point via private persons is for free. The owner of old appliances is responsible for bringing the appliances to these collecting points or to similar collection points. With this little personal effort, you contribute to recycle valuable raw materials and the treatment of toxic substances.

Caution:

- If you experience skin irritation when using the body massager, please refrain from using it, and wait two to three days to see if your symptoms ease. If symptoms persist or worsen, please consult a physician.
- If you have eczema, allergies or asthma, you may be more likely to experience a skin irritation or allergy from using the massager.
- The appliance it must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
- * This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- * Children shall not play with the appliance.
- * Cleaning and user maintenance shall not be made by children without supervision.

材质: 120g双铜+双面印刷

尺寸: 90*130mm

Pantone 3252C