
Instruction manual	Upute za uporabu
Bedienungsanleitung	εγχειρίδιο οδηγιών
Instrukcja obsługi	инструкция за употреба
Návod k použití	Manual de instructiuni
Návod na použitie	Instrukcijų vadovas
Használati utasítás	Kasutusjuhendit
Navodila	

Enjoy it.

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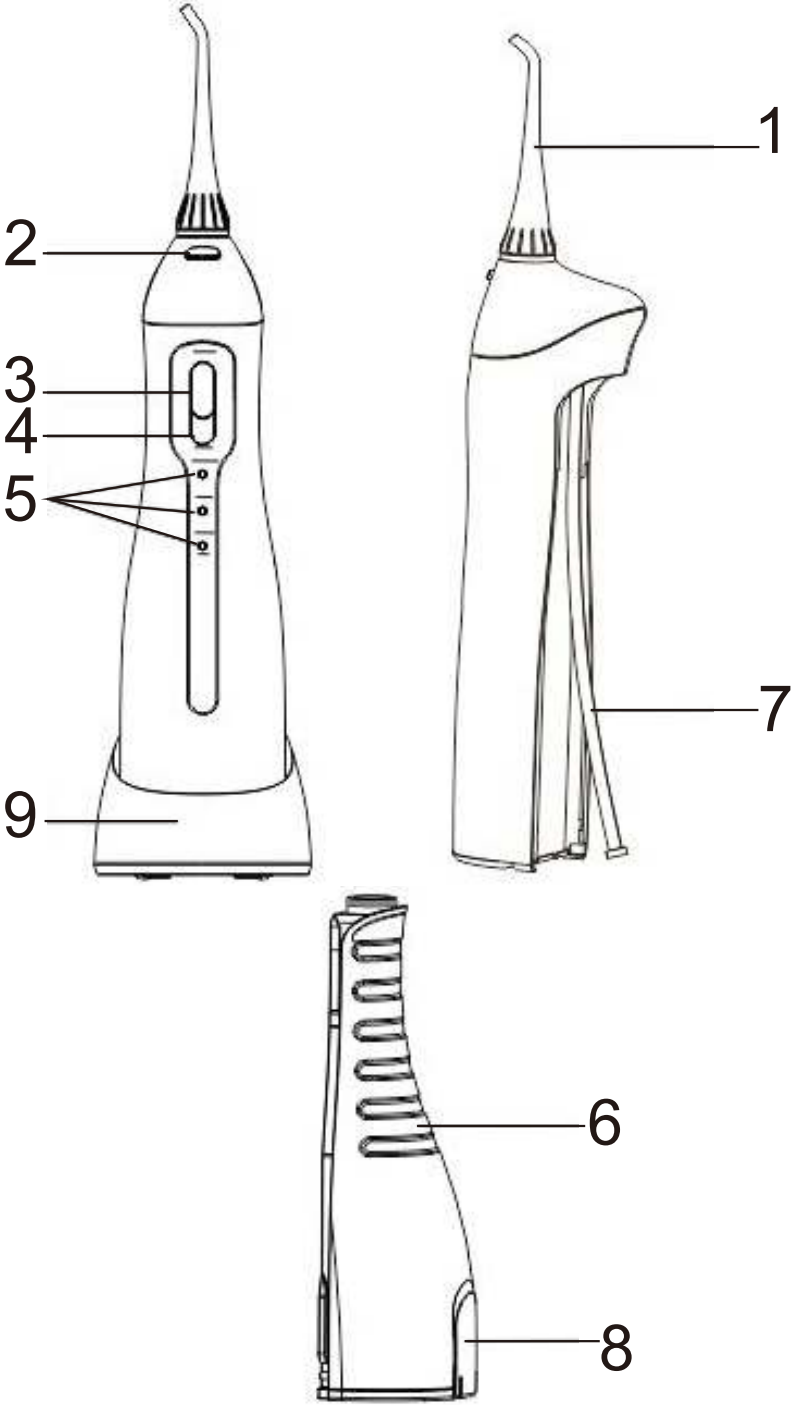Important notes

- Before operating the appliance, please read the instructions manual and follow the instructions provided.
- The manufacturer is not liable for any damages resulting from use of the appliance contrary to its intended use, or its inappropriate use. Save this instructions manual so you will also be able to refer to it when using the product in the future.
- The appliance is intended for home use only. Do not use for other purposes contrary to its intended use.
- Connect the appliance only to an electrical supply with specifications in accordance with the information on the data plate.
- This appliance may be used by children aged 8 years and above, and by persons with reduced physical and/or mental capabilities, and by persons with little experience and knowledge concerning the appliance, only if they have been given supervision or instructions concerning the use of the appliance in a safe way, and they understand the hazards involved. Children should not play with the appliance.
- Unattended children should not be allowed to perform the cleaning and maintenance of the appliance.
- Periodically check the state of the power cord. If the integral power cord is damaged, it should be replaced by the manufacturer, or a professional repair facility, or by a qualified person, in order to avoid hazard.
- The charging unit should be kept away from tanks containing water, such as bathtubs, washbasins or toilets.
- Use only the manufacturer's original spare parts, or those recommended by the manufacturer. Using accessories not recommended by the manufacturer can damage the

appliance and threaten operational safety.

- Do not pour hot water (hotter than 40 degrees Celsius), or any other corrosive liquid, into the container. It is not recommended to use mouthwash liquids with this irrigator, as some of them can damage the device.

CE Device is compliant with EU directives:

- Low voltage directive (LVD)
- Electromagnetic compatibility (EMC)

Device marked CE mark on rating label

When this crossed-out wheeled bin symbol is attached to a product it means that the product is covered by the European Directive 2012/19/EU. Please inform yourself about the local separate collection system for electrical and electronic products. Used appliance should be delivered to the dedicated collecting points due to hazardous components, which may affect the environment. Do not dispose this appliance in the common waste bin.

Your product contains batteries covered by the European Directive 2013/56/UE, which cannot be disposed with normal household waste. Please inform yourself about the local rules on separate collection of batteries because correct disposal helps to prevent negative consequences for the environmental and human health.

Environmental information

All unnecessary packaging has been omitted. We have tried to make the packaging easy to separate into three materials: cardboard (box), polystyrene foam (buffer) and poly ethylene (bags, protective foam sheet). Your system consists of material which can be recycled and reused if disassembled by a specialized company. Please observe the local regulations regarding the disposal of packaging materials, exhausted batteries and old equipment.

Your products is designed a manufactured with high quality materials and components, which can be recycled and reused.

Technical data

Power supply: 100-240V~50/60Hz

Power: 0.7-1.7W

- | | |
|-------------------------------------|----------------------|
| 1. Replaceable irrigation nozzles | 6. Water container |
| 2. Irrigation-nozzle release button | 7. Water-supply hose |
| 3. Start button | 8. Water-fill cap |
| 4. Operating-mode change button | 9. Charging station |
| 5. Operating-mode indicator | |

BATTERY CHARGING

You should make sure the irrigator is turned off before you place it in the charging station. Plug the charging station (9) into a power source. Place the irrigator in the charging station. The full charging process can take up to 24 hours. It is recommended to keep the irrigator in the charging station. The battery includes an overcharging-protection system. When the device is charging its battery, the charging indicator blinks. When the battery is fully charged, the indicator shines changelessly.

USAGE

1. Put the replaceable nozzle (1) onto the body.
2. Fill the water container. The water amount can be supplemented in two ways: by moving the inlet cap (8), or by moving the whole water container down. Make sure the water temperature is not higher than 40 °C
3. Set the mode with the operating-mode change button
4. Keep the device in the vertical position and direct the nozzle onto the teeth and inter-dental spaces.
5. Keep your mouth open so that excess water can go out.
6. Push the start button (3) to clean the teeth with water.
7. When cleaning is over, turn off the device with the start button (3). It is not recommended to leave water in the container, due to bacteria and other microorganism development.
8. To remove the nozzle, push the nozzle-release button (2) and pull the nozzle up.
9. Place the device's body in the charging station.

CLEANING MODES

1. NORMAL - standard cleaning mode
2. SOFT - delicate cleaning mode
3. PULSE - pulsing mode for gum stimulation

CLEANING AND MAINTANANCE

Before cleaning the device, make sure that the charging station is disconnected from the power supply.

Wipe the charging station with a damp cloth and dry it.

Wipe the irrigator's body with a damp cloth and dry it.

Do not submerge the irrigator's station nor its body in water.

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- CE** Device is compliant with EU directives:
- Low voltage directive (LVD)
 - Electromagnetic compatibility (EMC)
- Device marked CE mark on rating label

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All rights reserved. All brand names are registered trademarks of their